

2019  
NEW

# NAVIGATOR

ROLLATOR WITH FOREARM SUPPORTS



## ROTATING FOREARM SUPPORTS

The unique Navigator is a rollator with forearm supports specially designed for users with arthritis and/or limited strength in their arms and hands. It is equipped with very stable and comfortable height adjustable forearm supports. From now on, the Navigator is also available as a new version with horizontal angle-adjustable forearm supports that provides optimal adjustment options, so each user can set it exactly according to his needs and then enjoy the comfortable hand position. For users who need a rollator with forearm supports, Navigator is the best choice. For all the details and benefits, please read the following pages.



REHASENSE

# COMMON DETAILS & SPECIFICATIONS



The height adjustable forearm supports are very easy to adjust in a range of 75 – 115 cm (total for both rollator sizes).




Centre-to-centre distance between forearm supports is 47 cm. It is equipped with a strap in order to stabilize and strengthen the upper part. The strap must be connected during use.



The deep and soft padded forearm supports provides optimal arm comfort and support for the user. They are made of PU and the length is 24 cm. The surface is washable and easy to keep clean.



The brake and push handles are fully integrated into the forearm support and length adjustable in a range of 13 cm.

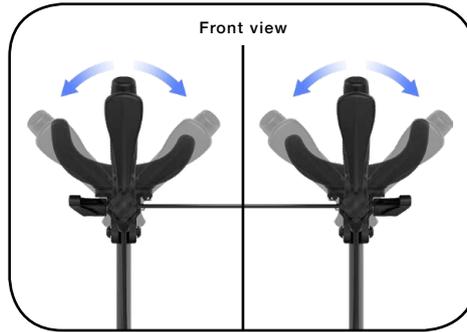


The push handles can be rotated to fit the most comfortable position of the user's hand. Width between push handles is 47 cm.



The hand brakes include parking brake. They run very smooth and can be operated with very little force. Reliable and maintenance free.



Easy to lock and unlock the rollator with the unique A+ locking system. It folds easily by pulling the strap on the seat. It stands firmly on its wheels taking up little space.

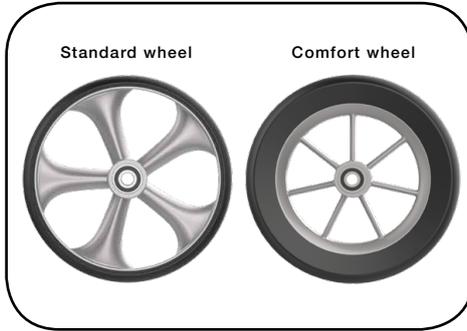

The Navigator is available with standard or comfort wheels. Standard wheels are made with black rubber tyre (TPE) and comfort wheels are made with black soft PU tyre.




Practical net bag as an extra feature to put between the handgrips for carrying small items such as a mobile phone, medicines, glasses, keys etc.



The Navigator is equipped with a step pedal on both sides which makes it easier to pass doorsteps and similar obstacles.




The specially designed aluminium frame is well-known of the other Rehasense rollator designs. The Navigator comes in 3 colours. For more information please turn to next pages.

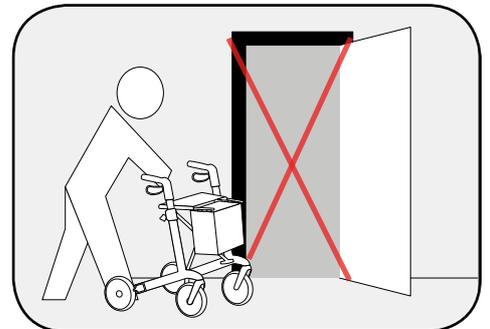

The Navigator is designed and approved for use on flat surfaces indoor only.

# NAVIGATOR WITH ROTATING FOREARM SUPPORTS

Provides optimal comfort & unique adjustment options (recommended)

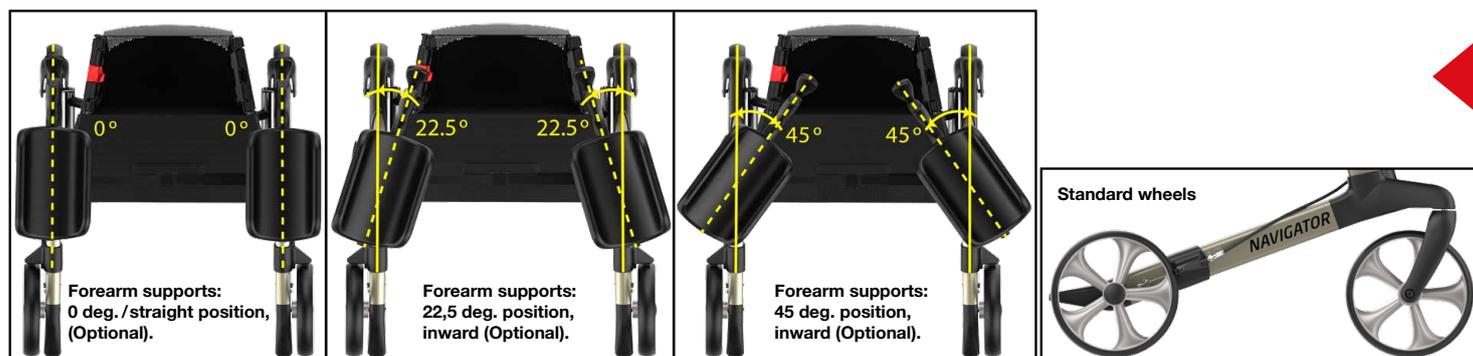


## NAVIGATOR WITH STANDARD WHEELS - Ø200 x 35 mm with black rubber tyre (TPE)

Color	Size / seat height	Seat width	Total weight w/o bag	Item number	Max user weight	Total width	Total length	Folded width	Push handle height	Height adjustment Dif. G/S
Champagne	L/62 cm	46 cm	9,9 kg	NRRCHL60	150 kg	63 cm	67 cm	28 cm	90-114 cm	10 steps

## NAVIGATOR WITH COMFORT WHEELS - Ø200 x 35 mm with black, high and soft PU tyre

Color	Size / seat height	Seat width	Total weight w/o bag	Item number	Max user weight	Total width	Total length	Folded width	Push handle height	Height adjustment Dif. G/S
Champagne	L/62 cm	46 cm	10,1 kg	NRRCHL60UP	150 kg	63 cm	67 cm	28 cm	90-114 cm	10 steps

# NAVIGATOR WITH FIXED FOREARM SUPPORTS

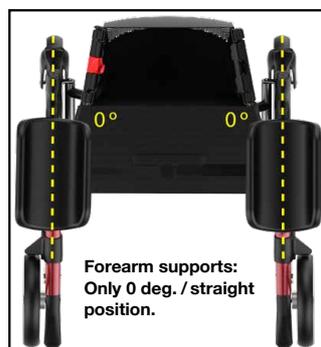

## NAVIGATOR WITH STANDARD WHEELS - Ø200 x 35 mm with black rubber tyre (TPE)

Color	Size / seat height	Seat width	Total weight w/o bag	Item number	Max user weight	Total width	Total length	Folded width	Push handle height	Height adjustment Dif. G/S
Sober grey	L/62 cm	46 cm	9,4 kg	NRSL600	150 kg	63 cm	67 cm	28 cm	90-115 cm	10 steps
Wine red	L/62 cm	46 cm	9,4 kg	NRWRL600	150 kg	63 cm	67 cm	28 cm	90-115 cm	10 steps
Sober grey	M/55 cm	46 cm	9,2 kg	NRSGM550	150 kg	63 cm	67 cm	28 cm	75-98 cm	7 steps

## NAVIGATOR WITH COMFORT WHEELS - Ø200 x 35 mm with black, high and soft PU tyre

Color	Size / seat height	Seat width	Total weight w/o bag	Item number	Max user weight	Total width	Total length	Folded width	Push handle height	Height adjustment Dif. G/S
Sober grey	L/62 cm	46 cm	9,6 kg	NRSL600UP	150 kg	63 cm	67 cm	28 cm	90-115 cm	10 steps
Wine red	L/62 cm	46 cm	9,6 kg	NRWRL600UP	150 kg	63 cm	67 cm	28 cm	90-115 cm	10 steps
Sober grey	M/55 cm	46 cm	9,4 kg	NRSGM550UP	150 kg	63 cm	67 cm	28 cm	75-98 cm	7 steps



## AVAILABLE COLORS:

WITH ROTATING FOREARM SUPPORTS



Champagne

NAVIGATOR WITH FIXED FOREARM SUPPORTS



Sober Grey



Wine Red

## ACCESSORIES

Rehasense rollator accessories are designed to make the use of your rollator more comfortable, practical, functional and safe during daily use. We are constantly focusing on innovation and development of new products that you can use to stay active. We value our customers and your feedback and ideas are always welcome. You can customize your rollator exactly to your need with a wide range of accessories as shown below. Ask your local dealer for a demonstration.




Net bag with zipper

Underseat net bag

Net bags




Crutch holder

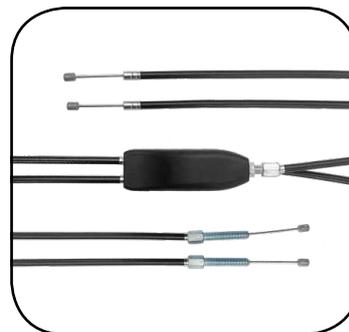

One-hand brake cable



Tray

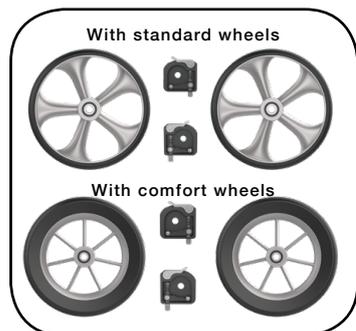

With standard wheels

With comfort wheels

Slow down brake



With zipper

With cover


Bags with lid



Infusion rack



Oxygen bottle holder

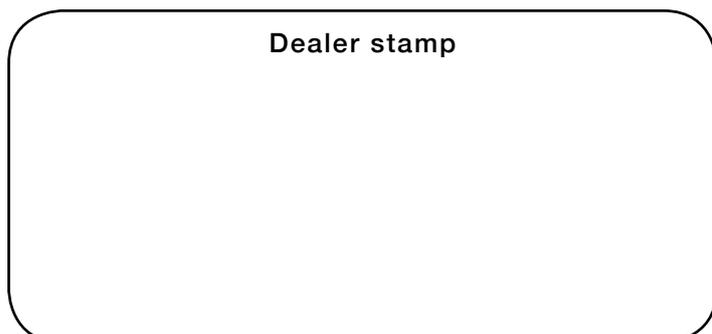

Dealer stamp



Tested & approved according to ISO 11199-3:2008

Ver 1 - 01 - 2019